

Certificate of Analysis

Nexus Research Laboratories EOOD

Nexus Research Laboratories EOOD
Str. Vrabcha No. 17A, Floor 1, Apt. 3,
Oborishte District
1000 Sofia
Bulgaria

Liquilabs s.r.o.

Inovační 122
252 41 Zlatníky-Hodkovice
Czechia
www.liquilabs.cz



[Verify Results Online](#)

Sample Identification

Sample Name Tesamorelin 10 mg
Batch Number XD-TSM10-PO05
Date Published 2026-09-19 09:47

Results for LYO-0887

Peptides	Result	Unit	Uncertainty	Acceptable Range
Tesamorelin Assay Peptide Screening 0.1% TFA	8.78	mg	[± 0.04]	
Tesamorelin Purity Peptide Screening 0.1% TFA	99.8	%	[± 0.5]	
Tesamorelin Identification by Spectrum Peptide Screening 0.1% TFA	981		[± 5]	
Tesamorelin Identification by RT Peptide Screening 0.1% TFA	0.999		[± 0.005]	
Mass Spectrometry	Result	Unit	Uncertainty	Acceptable Range
Molecular Ion Mass Identification (MS Deconvolution) Mass Spectrometry Identity	5135	Da	[± 1]	

	Method Specification		
Determination of identity, content and purity of Tesamorelin			
<i>Document number</i> TESA_008_2026	<i>Superseded document</i> -	<i>Number of pages</i> 4	

1. Content Assessment

1.1. Instrumentation

Module	Name	Serial Number
System Controller	Shimadzu CBM-40 Lite	L221226351398
Degassing Unit	Shimadzu DGU-403	NA
Pump	Shimadzu LC-40B XR	L22146350580
Autosampler	Shimadzu SIL-40C XR	L22216351622
Colum Thermostat	Shimadzu CTO-40S	L22236351602
PDA Detector	Shimadzu SPD-M40	L22276352808
SQ MS Detector	Shimadzu LCMS-2050	O12476200760

1.2. Chromatographic conditions

Chromatographic conditions	
Eluent A	0.1% DFA in Water (HPLC, Gradient Grade)
Eluent B	0.1% DFA in Acetonitrile (HPLC, Gradient Grade)
Flow rate	0.8 mL/min
Program	Gradient elution
Injection volume	2 µL
Colum Temperature	60°C
Column	Waters XSelect CSH C18, 100x2.1mm 2.5µm
Detection wavelength	280nm


Gradient Program		
Time [min]	A [%]	B [%]
1	98	2
13.99	35	65
14	5	95
15	5	95
15.01	98	2
19	end	

1.3. Sample preparation

Whole amount of container was dissolved in 2mL of water (LCMS Grade). 100 µL of sample was transferred to HPLC vial and diluted by 900 µL water (LCMS Grade) and submitted for analysis.

1.4. Calibration curve

Calibration curve detail	
Quantitative method	External Standard
Calibration Type	Linear
Number of calibration points	5
Force through Zero	Enabled
Weighting Method	None

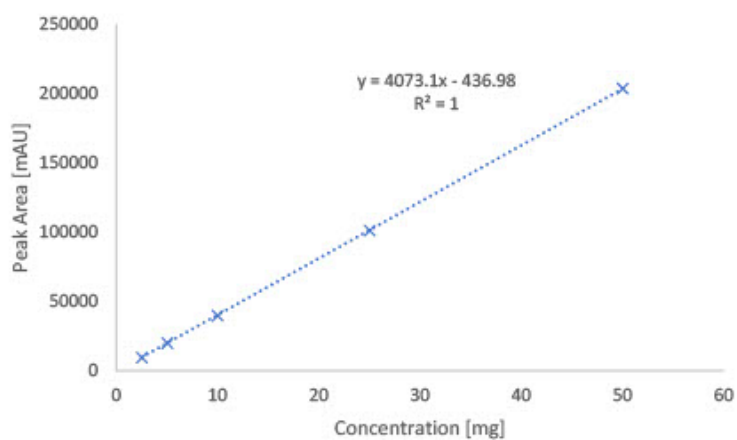

2. Purity assessment

2.1 Instrumentation

Module	Name	Serial Number
System Controller	Shimadzu CBM-40 Lite	L221226351398
Degassing Unit	Shimadzu DGU-403	NA
Pump	Shimadzu LC-40B XR	L22146350580
Autosampler	Shimadzu SIL-40C XR	L22216351622
Colum Thermostat	Shimadzu CTO-40S	L22236351602
PDA Detector	Shimadzu SPD-M40	L22276352808
SQ MS Detector	Shimadzu LCMS-2050	O12476200760

2.2 Chromatographic conditions

Chromatographic conditions	
Eluent A	0.1% DFA in Water (HPLC, Gradient Grade)
Eluent B	0.1% DFA in Acetonitrile (HPLC, Gradient Grade)
Flow rate	0.8 mL/min
Program	Gradient elution
Injection volume	2 µL
Colum Temperature	60°C
Column	Waters XSelect CSH C18, 100x2.1mm 2.5µm
Detection wavelength	225nm

Gradient Program		
Time [min]	A [%]	B [%]
1	98	2
13.99	35	65
14	5	95
15	5	95
15.01	98	2
19	end	

2.3 Purity assesment

Purity of compound assesed by area normalization method, comparing area of each peak to sum of area of all peaks detected at wavelenght of 225 nm.

3. Identity Assessment

3.1 Instrumentation

Module	Name	Serial Number
System Controller	Shimadzu CBM-40 Lite	L221226351398
Degassing Unit	Shimadzu DGU-403	NA
Pump	Shimadzu LC-40B XR	L22146350580
Autosampler	Shimadzu SIL-40C XR	L22216351622
Column Thermostat	Shimadzu CTO-40S	L22236351602
PDA Detector	Shimadzu SPD-M40	L22276352808
SQ MS Detector	Shimadzu LCMS-2050	O12476200760

3.2 Chromatographic conditions

Chromatographic conditions	
Eluent A	0.1% DFA in Water (HPLC, Gradient Grade)
Eluent B	0.1% DFA in Acetonitrile (HPLC, Gradient Grade)
Flow rate	0.8 mL/min
Program	Gradient elution
Injection volume	2 µL
Column Temperature	60°C
Column	Waters XSelect CSH C18, 100x2.1mm 2.5µm
Detection wavelength	280nm

Gradient Program		
Time [min]	A [%]	B [%]
1	98	2
13.99	35	65
14	5	95
15	5	95
15.01	98	2
19	end	

3.3 Molecular Ion Mass evaluation

Molecular ion mass was determined by deconvolution of multiply charged ESI-MS spectra to calculate the average neutral (zero-charge) molecular mass by equation:

$$M(\text{neutral}) = (z_i((mz_i) - H)) - ME$$

Where:

mz_i - Measured mass of charged particle

z_i - charge

H - proton mass (1.0076 Da)

ME - mass error

Analysis Report





Analysis of quantity of active ingredient by UHPLC with MS detection

Sample Information

Injection Volume : 2
Data File : LYO-0887__007.lcd
Method File : Peptide Screening_DFA_Group B.lcm
Date Acquired : 9/14/2026 7:00:35 PM

Chromatogram

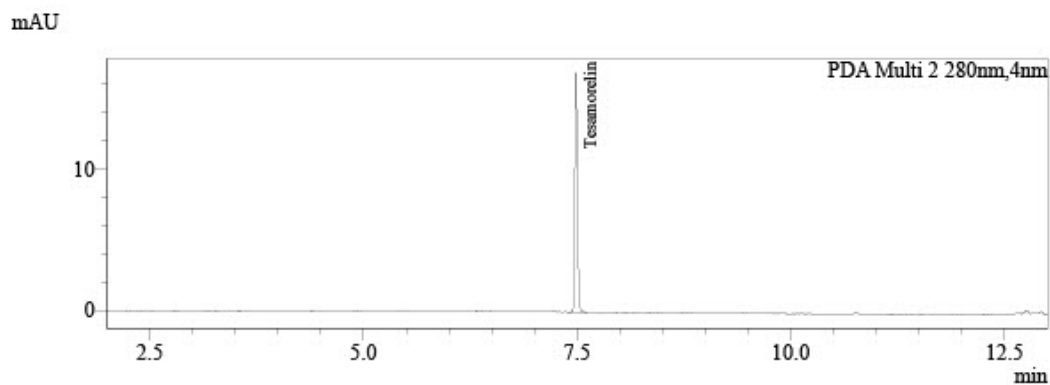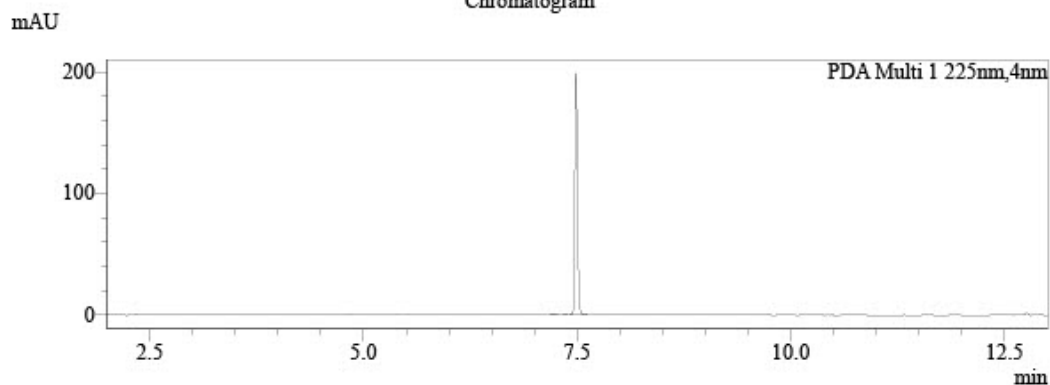

Peak Table

PDA Ch1 225nm

Peak#	Name	Ret. Time	Area	Area%	Conc.	Unit
1		7.220	768	0.183	0.000	
2		7.359	165	0.039	0.000	
3		7.485	418457	99.778	0.000	
Total			419389	100.000		

Peak Table


PDA Ch2 280nm

Peak#	Name	Ret. Time	Area	Area%	Conc.	Unit
1	Tesamorelin	7.485	35338	92.792	8.783	mg
2		14.220	2745	7.208	0.000	
Total			38083	100.000		

Analysis Report




Identification of active ingredient by UHPLC Mass Spectrometry

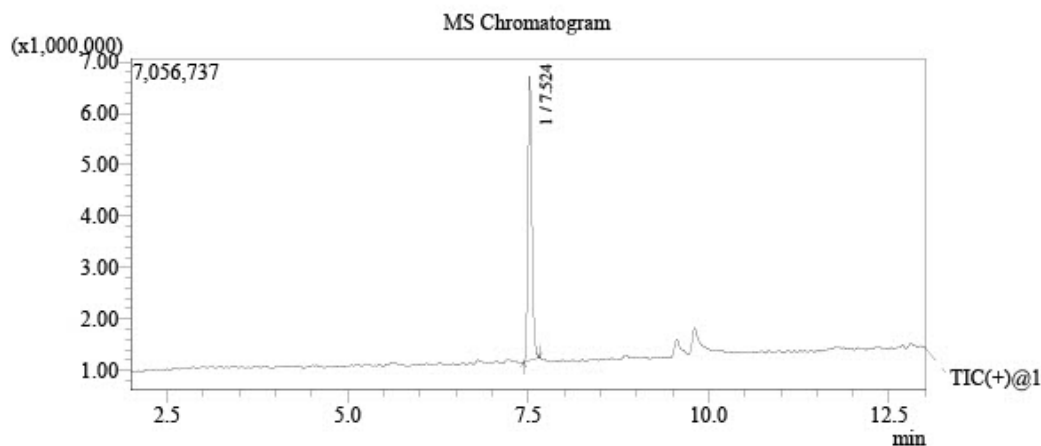


Library Search

Base Peak: 1047(942743)
Scan Range: 220.00 - 2000.00 Polarity: Positive

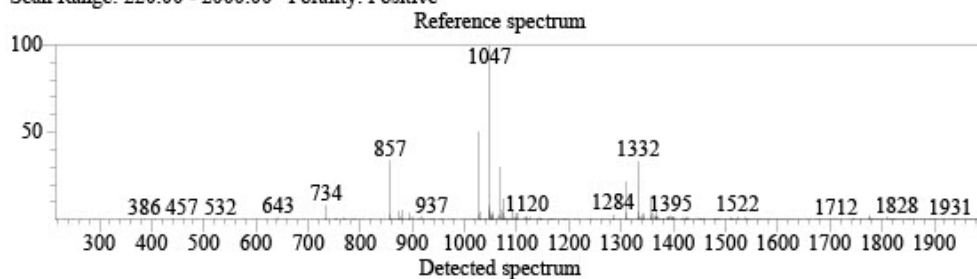

Library: Peptide Screening DFA.lib
Formula: C₂₂H₃₆N₇O₆S CAS: 218949-48-5 Mol.Weight: 5135(0!=)
Compound Name: Tesamorelin

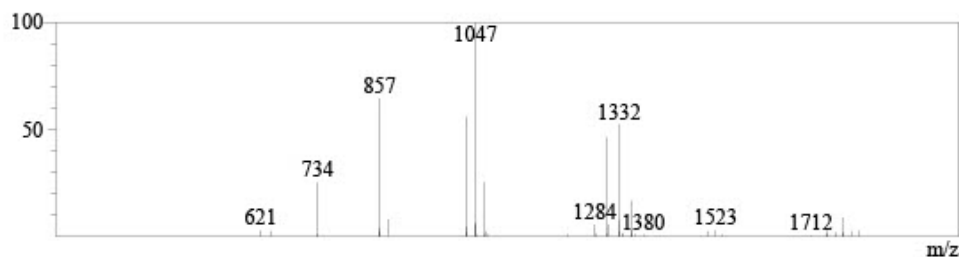

Responsibles



Mr. Ján Galbavý
CEO

Analysis results relate only to the samples tested.

This document shall not be reproduced except in full, without the written approval of Liquilabs s.r.o.